

Fast Determination of Inorganic Anions in Municipal Drinking Water Using Capillary Ion Chromatography

Fei Pang, Terri Christison, Richard Jack, and Linda Lopez, Sunnyvale, CA, USA

Introduction

The determination of common inorganic anions in drinking water is one of the most important applications of ion chromatography (IC) worldwide. IC has been approved for compliance monitoring of common inorganic anions in U.S. drinking water since the mid-1980s, as described in U.S. EPA Method 300.0.¹ A considerable number of regulatory IC methods have been published in many other industrialized countries (e.g., in Germany, France, Italy, and Japan). In addition, many standard organizations (including ISO, ASTM, and AWWA) have validated IC methods for the determination of inorganic anions in drinking water.^{2,3} The concentrations of some anions in drinking water are regulated due to their toxic effects. For example, high levels of fluoride cause skeletal and dental fluorosis, and nitrite and nitrate can cause methemoglobinemia, which can be fatal to infants. Other common anions, such as chloride and sulfate, are considered secondary contaminants and can affect odor, color, and certain aesthetic characteristics.

This application brief describes the fast determination of inorganic anions in drinking water samples with the Thermo Scientific Dionex ICS-5000 capillary system tolerating 5000 psi backpressure. This capability facilitates higher sample throughput by simply increasing the flow rate of analysis. Figure 1 shows the analysis of inorganic anions using a Thermo Scientific Dionex IonSwift MAX-100 capillary column and suppressed conductivity detection. The Dionex IonSwift™ MAX-100 column uses monolith technology and is designed for the separation of small molecules including organic acids and inorganic anions using a hydroxide gradient delivered by an eluent generator. All analytes of interests are baseline resolved and elute within 8 min at 24 $\mu\text{L}/\text{min}$ flow rate. At this flow rate, about 4500 psi backpressure is generated.

Scaling down from standard bore to capillary scale brings many benefits to IC users. One of the most important values is that the system can be always on and ready for analysis because of its low consumption of eluent (15 mL of source water a day at 10 $\mu\text{L}/\text{min}$ flow rate). The system remains stable when operated continuously, greatly reducing time spent equilibrating and calibrating. The amount of waste generated is significantly decreased and the Thermo Scientific Dionex EGC Eluent Generation Cartridge producing the eluent typically lasts 18 months under continuous operation, which translates into reduced overall cost of ownership.

Key Words

- Capillary IC
- HPIC
- Dionex ICS-5000
- U.S. EPA Method 300.0
- Dionex IonSwift MAX-100

Column:	Dionex IonSwift MAX-100G/ MAX-100, 0.25 mm	Samples: A: Bottled water B: UK Municipal Drinking Water C: Standard
Eluent:	Dionex EGC-KOH (Capillary)	
Gradient:	0.1 mM KOH (0.1 min), 0.2–2 mM (0.1 to 1 min), 2–15 mM (1 to 4 min), 15–35 mM (4 to 6 min), 35 mM (6 to 10 min)	Peaks:
Flow Rate:	24 $\mu\text{L}/\text{min}$	1. Fluoride 0.06 0.06 2 mg/L
Inj. Volume:	0.4 μL	2. Chloride 30.9 8.6 10
Temperature:	35 $^{\circ}\text{C}$	3. Nitrite — — 10
Detection:	Suppressed Conductivity, Thermo Scientific Dionex ACES 300 Anion Capillary Electrolytic Suppressor	4. Bromide — — 10
		5. Nitrate 12 3.5 10
		6. Carbonate — — 10
		7. Sulfate 104 13.5 10
		8. Phosphate 3.1 — 20

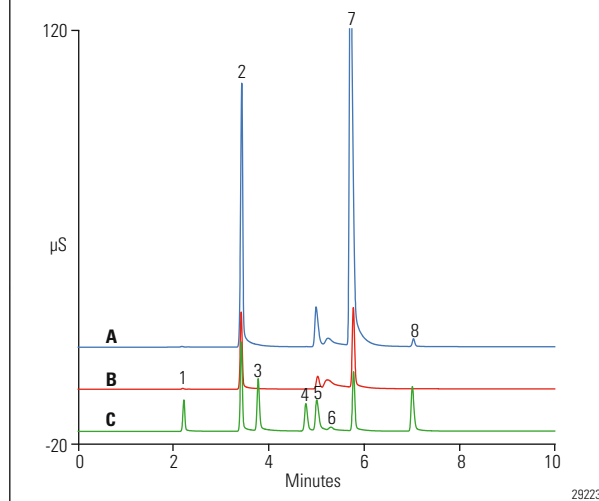


Figure 1. Fast determination of anions in two drinking water samples using high-pressure capillary IC.

29223

Conditions

Dionex ICS-5000 capillary system, Thermo Scientific Dionex AS-AP Autosampler, and Thermo Scientific Dionex Chromeleon Chromatography Data System software are used in this experiment. All experimental parameters are listed in Figure 1.

Sample Information and Preparation

Municipal drinking water was analyzed by capillary IC without sample pretreatment.

References

1. The Determination of Inorganic Anions in Water by Ion Chromatography; Method 300.0, Revision 2.1; U.S. Environmental Protection Agency; Cincinnati, Ohio, 1993.
2. Greenberg, A. E.; Clesceri, L. S.; Eaton, A. D., Eds.; Standard Methods for the Examination of Water and Wastewater, 18th ed.; Am. Public Health Assoc.; Washington, DC, 1992.
3. Standard Test Methods for Anions in Water by Chemically Suppressed Ion Chromatography; D4327-97, Vol. 11.01; American Society for Testing and Materials; West Conshohocken, Pennsylvania, 1999; p. 420-427.

In addition to these offices, Thermo Fisher Scientific maintains a network of representative organizations throughout the world.

Dionex Products:
(408) 737-0700

**North America:
U.S./Canada**
(847) 295-7500

**South America:
Brazil**
(55) 11 3731 5140

**Europe:
Austria**
(43) 616 51 25

Benelux
(31) 20 683 9768
(32) 3 353 4294

Denmark
(45) 36 36 90 90

France
(33) 1 39 30 01 10

Germany
(49) 6126 991 0

Ireland
(353) 1 644 0064

Italy
(39) 02 51 62 1267

Sweden
(46) 8 473 3380

Switzerland
(41) 62 205 9966

United Kingdom
(44) 1276 691722

**Asia Pacific:
Australia**
(61) 2 9420 5233

China
(852) 2428 3282

India
(91) 22 2764 2735

Japan
(81) 6 6885 1213

Korea
(82) 2 2653 2580

Singapore
(65) 6289 1190

Taiwan
(886) 2 8751 6655

www.thermoscientific.com

Legal Notices: ©2012 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. This information is presented as an example of the capabilities of Thermo Fisher Scientific Inc. products. It is not intended to encourage use of these products in any manners that might infringe the intellectual property rights of others. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.



Thermo Scientific Dionex products are designed, developed, and manufactured under an ISO 9001 Quality System.

AB141_E 02/12S LPN3059